



# BJOY Hand- C

## Access to computer with minimal finger strength



Ref:  
BJ-884-C



The BJOY Hand-C allows you to control your computer with the slightest movement of your fingers. Designed to be comfortable and ergonomically used with one hand, the BJOY Hand-C is based on a very sensitive microswitch joystick that requires minimal strength. The BJOY Hand C is very similar to the BJOY Hand A but its microswitch joystick requires less accuracy to be used. It is part of the BJOY range of products, which means its hardware can be tailored to each user's needs and capabilities and its operation customized using the BJOY Configurator Software.


- ✓ Use your mouse **with precision** with **minimal** finger strength.
- ✓ Get to your favorites applications in an **easy and simple way**.
- ✓ **Personalize** your mouse so it is more suited to your needs. We will make it at your request, based on your individual needs
- ✓ **"Plug and Play"** devices, compatible with Windows, Mac OS and Linux.
- ✓ Use the **BJOY Configurator** software to modify each button's function, speed, orientation and many other features.
- ✓ All settings adjusted with the BJOY Configurator are **stored in the device** so you can take them with you wherever you go.



### TECHNICAL FEATURES

- "Plug and Play" USB device.
- Size: 50 x 75 x 17 mm .
- Four small ultrasensitive push buttons with customizable functions.
- Microswitch joystick: movements in 8 directions.
- Speed independent from the action performed on the joystick.
- Very sensitive joystick.
- Ergonomically designed for one-hand operation.

### CUSTOMIZE THE OPERATION SETTINGS OF YOUR MOUSE WITH THE BJOY CONFIGURATOR!

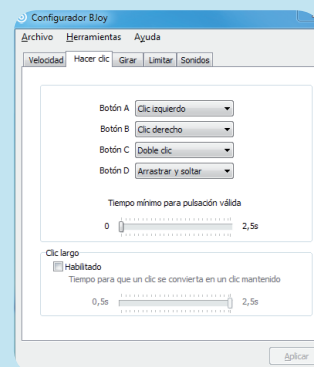

#### THE BJOY CONFIGURATION SOFTWARE ALLOWS YOU TO:

- Adjust the speed, sensitivity, acceleration and auditory feedback of the joystick according to the user's needs.
- Orientate the device easily so you can adjust it to the user's position.
- Each button can perform the following actions: left click, right click, drag & drop, double click or nothing.
- Determine the time needed for a button press to be considered valid.
- Configure the Drag & drop function as a long press on the button you select.
- Configure the joystick with 4 or 8 directions to simplify access for people with inaccurate movements.

**BJOY ALLOWS YOU TO CREATE YOUR OWN  
CUSTOM MADE MOUSE!**